



MS-H630

Auto Hematology Analyzer

Accessible Utilities In Optimizing Hematology





Accurate and reliable



Bi-directional interface



User-friendly




Safe and secure



User-friendly interface (10.4"Intuitive color touch screen menus; simple operation, virtual keyboard)



Lyse reagent internal placement, non-cyanide method

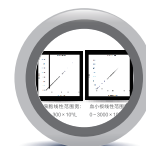

Excellent reproducibility and broad linearity range



Integrated thermal printer with comprehensive data management, compatible with USB external printer



Product specification

Throughput	Sample Volume		Storage
60 Samples/ hour	CBC 9μL		50,000 Sample results, including parameter, histogram, and patient information
Analysis Channel	Counting modes		Input/output
Dual-Channel for counting	Venous whole blood, capillary whole blood, and pre-diluted blood		LAN connection and 4 USB ports for mouse, keyboard, barcode scanner
Innovative Technology	Reference Value Setting		Operating environment
Automatic floating discriminators, Plus baseline tracking technology, Waveform discrimination technology, Digital Sheath-Flow technology	Pediatrics, Male/Female adult, general and user-defined		Temperature: 10°C ~ 30°C Humidity: 20% ~ 85% Air pressure: 70kPa~106kPa
Principle	Flagging system		Power requirement
Impedance method (WBC/RBC/PLT) and cyanide free colorimetric method (HGB)	Reliable alarm system with customized function		A.C100V~240V (50/60Hz), ≤200VA
Parameters	Result Performance		Carryover
21 Reportable Parameters CBC: Whole Blood Mode: WBC, Lym#, Mid#, Gran#, Lym%, Mid%, Gran%, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV, PDW, PCT, P-LCR, P-LCC 3 Histograms: WBC (3-Part Differential), RBC, PLT Sensitive flagging (WBC, RBC & PLT) for Pathological sample and counting failure flags	Linearity	Precision(CV%)	WBC, RBC, HGB ≤0.5%, PLT ≤1.0%
	WBC 0~300 x 10 ⁹ /L	≤2.0%(7.0~15.0) ≤3.5%(4.0~6.9)	Dimensions
	RBC 0~8.5 x 10 ¹² /L	≤1.5%(3.5~6.5)	
	HGB 0~250g/L	≤1.5%(100~180)	Weight
	HCT 0~67%		
	PLT 0~3000 x 10 ⁹ /L	≤4.0%(150~500)	

MEDICALSYSTEM BIOTECHNOLOGY CO., LTD.

ADD: No.299, Qiming South Road, Yinzhou District, Ningbo City, Zhejiang Prov. China
 TEL: +86-574-55337009 FAX: +86-574-88231889
 ZIP Code: 315104 E-mail: sales@nbmedicalsystem.com



www.nbmedicalsystem.com